

Acura® manual
your expectations – our tradition

new
Qualitix®
tip line

digital reading micropipettes

Acura® manual pipettes are precision instruments combining high-tech materials with exceptional, user friendly ergonomics. A wide choice of models offers much more than superior performance and makes pipetting safe and comfortable.

The Acura® manual line features:

- Effortless activation
- Smart and reliable volume adjustment
- Precision digital display visible at all time
- *Justip*™ adjustable tip ejector fitting most tips
- *Swift-set* user calibration system
- Shock, UV-light and autoclaving resistance



The Acura® manual models

825 **826** **XS**
micro micro

815 **835** **855**
fixed macro multi

 **SOCOREX**
SWISS



► Most comfortable pipetting

The line fits everybody's hand with or without glove. Ergonomic shape, light weight and soft plunger travel care for optimal comfort.



► Smart and reliable volume adjustment

State-of-the-art micrometric setting is performed smoothly and precisely by turning the push button. Finely machined click-stops and free rotating cap prevent any unwanted volume alteration.

► Digital display always visible at a glance

Located on the front of the instrument, the display window allows the user to comfortably read the set volumes at any time during pipetting. A comforting reassurance for the user.

Universal micropipettes

The classic Acura® manual 825 is intended for all applications. It brings convenience, comfort, safety and sturdiness. A wise choice for your laboratory.

825
micro

0.1 – 2 µL
0.5 – 10 µL
1 – 10 µL
2 – 20 µL
5 – 50 µL
10 – 100 µL
20 – 200 µL
100 – 1000 µL



Work station Twister™ universal 336

Rotary stand holds up to six single channel pipettes of most brands. Seven translucent colours to choose from.



► Acura® manual 825/835 Triopack

Nine different packs up to 10 mL, each containing three adjustable micropipettes, to cover all individual requirements in any research or routine application. Available at attractive budget prices, they match with today's budget constraints.

In-line™ universal 337

Space saving stand holds up to seven single channel pipettes of most brands. Anti-slip pads guarantee bench stability. Available in five pastel colours.

Extra sharp micropipettes

Derived from the Acura® *manual* 825, the 826 XS line was designed with the research scientists in mind. Shorter, smoother, lighter in weight, the eight models offer key features aiming to excellence in pipetting.

826 XS micro extra sharp

0.1 – 2 µL
0.5 – 10 µL
1 – 10 µL
2 – 20 µL
5 – 50 µL
10 – 100 µL
20 – 200 µL
100 – 1000 µL



Short shaft



Narrow conical shaft end

Small diameter offers easy access to microtubes.

Short shaft

Length reduction greatly enhances precision of hand movement and guarantees superior drivability.

A must for research scientists

The extra sharp pipette is optimal for any fine applications such as pipetting in microtubes or microplates.

Performance and weight Acura® *manual* 825 and 826 XS

Volume	Weight		Inaccuracy (E%)			Imprecision (CV%)		
	825	826	Min. vol.	Mid. vol.	Max. vol.	Min. vol.	Mid. vol.	Max. vol.
0.1 – 2 µL	88 gr	84 gr	<± 6.0 % ¹⁾	<± 4.0 %	<± 2.0 %	< 5.0 % ¹⁾	< 3.3 %	< 1.5 %
0.5 – 10 µL	88 gr	84 gr	<± 2.5 % ²⁾	<± 1.8 %	<± 1.0 %	< 1.8 % ²⁾	< 1.2 %	< 0.5 %
1 – 10 µL	88 gr	85 gr	<± 2.5 %	<± 1.8 %	<± 1.0 %	< 2.5 %	< 1.6 %	< 0.7 %
2 – 20 µL	88 gr	85 gr	<± 2.5 %	<± 1.8 %	<± 1.0 %	< 1.7 %	< 1.0 %	< 0.5 %
5 – 50 µL	90 gr	86 gr	<± 1.5 %	<± 1.3 %	<± 1.0 %	< 1.0 %	< 0.7 %	< 0.4 %
10 – 100 µL	92 gr	87 gr	<± 1.5 %	<± 1.2 %	<± 0.8 %	< 1.0 %	< 0.6 %	< 0.2 %
20 – 200 µL	90 gr	86 gr	<± 1.5 %	<± 1.1 %	<± 0.6 %	< 0.6 %	< 0.4 %	< 0.2 %
100 – 1000 µL	95 gr	88 gr	<± 1.5 %	<± 1.0 %	<± 0.5 %	< 0.5 %	< 0.4 %	< 0.2 %


Performance values obtained with bidest. water at constant temperature (± 0.5°C) comprised between 20 and 25°C in accordance with ISO 8655.

¹⁾ measured at 0.5 µL ²⁾ measured at 1 µL



► Tip ejection - easier than ever

Most ergonomically positioned, the large surface, soft padded ejector button ensures low-pressure activation. In addition to ease of tip ejection, the shaft height adjustment system called *Justip™* (4 mm span) - controlled by efficient click-stops - allows a wide selection of tips to tightly fit the nozzle.

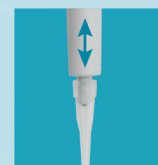

► Extra smooth activation

Unique tightness lip seal and new springs provide ultra-soft pipetting, reducing hand fatigue while working.

Pipetting		Overshoot	
Start	End	Start	End
≤1.6 N*	≤2.8 N*	≤9.3 N*	≤11.5 N*


* 1 Newton (N) ~ 0.1 kilogramme force (kgf)

Above forces, measured on a 20 - 200 µL model, are indicative of very limited finger efforts. Yet, the overshoot stop offers a clear tactile indicator.



► Acura® *manual* XS 826 TwixS pack

Six different packs up to 1000 µL, each containing two adjustable micropipettes and a free pipette holder, to cover all individual requirements in any research application. Available at attractive prices, they match with today's budget constraints.

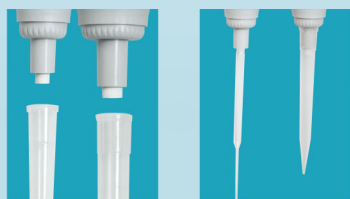

► Nozzle protection filter

Macro models accommodate interchangeable nozzle filters for added protection against liquid overflow and barrel contamination.

► Pasteur pipette adapter

Optional adapters allow fitting of glass Pasteur pipette on 2 and 5 mL models in addition to Polypropylene tips. The adapters include two O-rings guaranteeing proper holding and tightness of any standard 2 mL Pasteur pipettes (Ø 7 mm).



► Smarties caps colour coding at will

Free rotating colour caps protecting against unwanted volume alteration, can be changed for individual, lab or department identifications. Pick among 14 different colours or choose the smartie mix including all of them for unlimited combinations.



► Shelf pipette holder 332 for two instruments

Fitting micro, macro and/or multi-channel pipettes, its attractive design accommodates all Acura[®] and Calibra[®] instruments. The self-adhesive strip supplied ensures lasting attachment on any surface.

Fixed volume micro and macropipettes

They show most stable performance and provide consistent results for any analytical or routine tests.

No risk of erroneous volume setting.

Available in a selection of micro sizes from 1 to 1000 µL and macro volumes models of 2, 2.5, 5 and 10 mL.

815 / 835F

fixed

1 µL to 10 mL



Multi and single channel pipette work station 340

Innovative design fits three Calibra[®] and Acura[®] models. Easy reach of each instrument. Heavy bottom plate and anti-slip pads guarantee high stability.

Performance and weight Acura[®] manual 815 and 835F

Volume	Weight	Inaccuracy (E%)	Imprecision (CV%)	Volume	Weight	Inaccuracy (E%)	Imprecision (CV%)
Acura[®] manual 815							
1 µL	85 gr	<± 2.0 %	< 1.0 %	150 µL	89 gr	<± 0.7 %	< 0.3 %
5 µL	85 gr	<± 1.4 %	< 1.0 %	200 µL	89 gr	<± 0.7 %	< 0.3 %
10 µL	85 gr	<± 0.7 %	< 0.7 %	250 µL	93 gr	<± 0.7 %	< 0.4 %
10Y µL	85 gr	<± 0.7 %	< 0.8 %	300 µL	93 gr	<± 0.7 %	< 0.4 %
15 µL	85 gr	<± 0.7 %	< 0.5 %	400 µL	93 gr	<± 0.7 %	< 0.3 %
20 µL	85 gr	<± 0.7 %	< 0.5 %	500 µL	93 gr	<± 0.7 %	< 0.3 %
25 µL	89 gr	<± 0.7 %	< 0.6 %	600 µL	93 gr	<± 0.7 %	< 0.3 %
30 µL	89 gr	<± 0.7 %	< 0.6 %	700 µL	93 gr	<± 0.7 %	< 0.2 %
32 µL	89 gr	<± 0.7 %	< 0.6 %	750 µL	93 gr	<± 0.7 %	< 0.2 %
40 µL	89 gr	<± 0.7 %	< 0.5 %	800 µL	93 gr	<± 0.7 %	< 0.2 %
50 µL	89 gr	<± 0.7 %	< 0.4 %	900 µL	93 gr	<± 0.7 %	< 0.2 %
60 µL	89 gr	<± 0.7 %	< 0.5 %	1000 µL	93 gr	<± 0.6 %	< 0.2 %
70 µL	89 gr	<± 0.7 %	< 0.4 %	Acura[®] manual 835F			
75 µL	89 gr	<± 0.7 %	< 0.4 %	2 mL	105 gr	<± 0.5 %	< 0.2 %
80 µL	89 gr	<± 0.7 %	< 0.35 %	2.5 mL	111 gr	<± 0.8 %	< 0.5 %
90 µL	89 gr	<± 0.7 %	< 0.3 %	5 mL	111 gr	<± 0.6 %	< 0.3 %
100 µL	89 gr	<± 0.7 %	< 0.3 %	10 mL	120 gr	<± 0.5 %	< 0.2 %
120 µL	89 gr	<± 0.7 %	< 0.4 %				

Performance values obtained with bidest. water at constant temperature (± 0.5°C) comprised between 20 and 25°C in accordance with ISO 8655.

Macropipettes

Optimal for in- and outdoor environmental analysis, clinical chemistry and cell culture. They offer superior performance and make pipetting safer than ever.

835
macro

0.2 – 2 mL
0.5 – 5 mL
1 – 10 mL



Multichannel pipettes

Instruments extend pipetting possibilities in 96-well microplates and guarantee outstanding performance and results.



855
multi

0.5 – 10 µL
5 – 50 µL
10 – 100 µL
20 – 200 µL
40 – 350 µL




► Optimal working position

In addition to the ergonomic design and perfect hand fitting of the pipette, the multichannel volumetric module revolves to allow the selection of the best working position.



► Adjustable tip ejector

Height of tip ejector can be adjusted in no time within a 4 mm span to perfectly match to tip used. Bowed shape of ejector head guarantees effortless, easy sequential tip ejection.

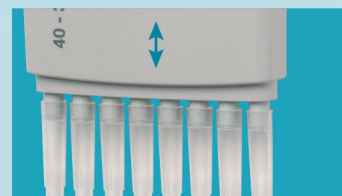

► All purpose reservoirs

Well adapted for multichannel pipettes, the reagent reservoirs offer various shapes in different configurations at economical cost. See page 10.

Performance and weight Acura® manual 835

Volume	Weight	Inaccuracy (E%)			Imprecision (CV%)		
		Min. vol.	Mid. vol.	Max. vol.	Min. vol.	Mid. vol.	Max. vol.
0.2 – 2 mL	108 gr	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %
0.5 – 5 mL	116 gr	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %
1 – 10 mL	130 gr	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %	<± 1.5 %

Performance and weight Acura® manual 855

Volume	Weight		Inaccuracy (E%)			Imprecision (CV%)		
	8x	12x	Min. vol.	Mid. vol.	Max. vol.	Min. vol.	Mid. vol.	Max. vol.
0.5 – 10 µL	155 gr	184 gr	<± 3.5 % ¹⁾	<± 3.5 % ¹⁾	<± 3.5 % ¹⁾	<± 3.5 % ¹⁾	<± 3.5 % ¹⁾	<± 3.5 % ¹⁾
5 – 50 µL	165 gr	199 gr	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %
10 – 100 µL	165 gr	165 gr	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %
20 – 200 µL	165 gr	201 gr	<± 0.9 %	<± 0.9 %	<± 0.9 %	<± 0.9 %	<± 0.9 %	<± 0.9 %
40 – 350 µL	168 gr	203 gr	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %	<± 1.0 %

Performance values obtained with bidest. water at constant temperature (± 0.5°C) comprised between 20 and 25°C in accordance with ISO 8655. ¹⁾ measured at 1 µL

Ordering information - Qualitix® pipette tip selection

10 µL Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
  	31.4	Natural	10 µL ultra-microtips, fine pointed, graduated, autoclavable		1000/bag	10 bags	3090.0010B
		Natural	10 µL ultra-microtips, fine pointed, graduated, autoclavable		12x 96/rack	10 packs	3090.0010R
		Natural	10 µL <i>tipfill™</i> system incl. ultra-microtips, fine pointed, graduated		13x 96/pack	10 packs	3090.0010T
			Empty rack for 10/200 µL tips to be used with <i>tipfill™</i> system, autoclavable		10/pack	10 packs	3070.0210ER
		Natural	10 µL ultra-microtips, fine pointed, graduated, autoclavable	✓	12x 96/rack	10 packs	3090.0010RS
	38.1	Natural	10 µL extended microtips, fine pointed, autoclavable		1000/bag	10 bags	3020.0010B
		Natural	10 µL extended microtips, fine pointed, autoclavable		10x 96/rack	10 packs	3020.0010R
	71.0	Natural	10 µL capillary tips for gel loading, autoclavable		4x 204/rack	N/A	302.0010GR

10 µL Filter Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	31.4	Natural	10 µL ultra-microtips, fine pointed, graduated	✓	12x 96/rack	10 packs	3090.0010FRS
	38.1	Natural	10 µL extended microtips, fine pointed	✓	10x 96/rack	10 packs	3020.0010FRS

10 µL Low Protein Binding Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	31.4	Natural	10 µL ultra-microtips, low binding, fine pointed, autoclavable		12x 96/rack	10 packs	3090.0010PR
		Natural	10 µL ultra-microtips, low binding, fine pointed, autoclavable	✓	12x 96/rack	10 packs	3090.0010PRS
	38.1	Natural	10 µL ultra-microtips, low binding, fine pointed, autoclavable		10x 96/rack	10 packs	3020.0010PR
		Natural	10 µL ultra-microtips, low binding, fine pointed, autoclavable	✓	10x 96/rack	10 packs	3020.0010PRS

20 - 100 µL Filter Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	49.6	Natural	20 µL filter tips, bevel point	✓	10x 96/rack	10 packs	3080.0020FRS
	49.6	Natural	100 µL filter tips, bevel point	✓	10x 96/rack	10 packs	3080.0100FRS

200 µL Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
  	49.6	Natural	200 µL universal tips, bevel point, graduated, autoclavable		1000/bag	10 bags	3070.0200B
		Natural	200 µL universal tips, bevel point, graduated, autoclavable		10x 96/rack	10 packs	3070.0200R
		Natural	200 µL <i>tipfill™</i> system incl. universal tips, bevel point, graduated		10x 96/pack	10 packs	3070.0200T
			Empty rack for 10/200 µL tips to be used with <i>tipfill™</i> system, autoclavable		10/pack	10 packs	3070.0210ER
		Natural	200 µL universal tips, bevel point, graduated, autoclavable	✓	10x 96/pack	10 packs	3070.0200RS
 	49.6	Yellow	200 µL regular tips, bevel point, autoclavable		1000/bag	10 bags	3290.0200B
		Yellow	200 µL regular tips, bevel point, autoclavable		10x 96/rack	10 packs	3290.0200R
		Yellow	200 µL regular tips, bevel point, autoclavable	✓	10x 96/rack	10 packs	3290.0200RS
 	49.0	Yellow	200 µL short neck tips for Calibra® pipettes, autoclavable		2x 500/bag	N/A	328.0200B
		Yellow	200 µL short neck tips for Calibra® pipettes, autoclavable		10x 96/rack	N/A	328.0200R
	82.5	Natural	200 µL capillary tips for gel loading, autoclavable		6x 96/rack	8 packs	3080.0200GR

200 µL Filter Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	51.2	Natural	200 µL universal tips, bevel point	✓	10x 96 / rack	10 packs	3070.0200FRS
	83.8	Natural	200 µL extended tips, fine pointed	✓	6x 96 / rack	8 packs	3080.0200LFRS

200 µL Low Protein Binding Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	49.6	Natural	200 µL universal tips, low binding, graduated, autoclavable		10x 96 / rack	10 packs	3070.0200PR
		Natural	200 µL universal tips, low binding, graduated, autoclavable	✓	10x 96 / rack	10 packs	3070.0200PRS

350 µL Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	58.8	Natural	350 µL tips, autoclavable		1000 / bag	N/A	308.0350B
		Natural	350 µL tips, autoclavable		10x 96 / rack	N/A	308.0350R

300 µL Filter Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	58.8	Natural	300 µL tips	✓	10x 96 / rack	N/A	308.0300FR

1000 µL Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	88.9	Natural	1000 µL universal tips, fine pointed, autoclavable		1000 / bag	10 bags	3070.1000B
		Natural	1000 µL universal tips, fine pointed, autoclavable		6x 96 / rack	8 packs	3070.1000R
		Natural	1000 µL <i>tipfill™</i> system incl. universal tips, fine pointed		5x 96 / rack	10 packs	3070.1000T
			Empty rack for 1000 µL tips to be used with <i>tipfill™</i> system, autoclavable		6 / pack	8 packs	3070.1000ER
		Natural	1000 µL universal tips, fine pointed, autoclavable	✓	6x 96 / rack	8 packs	3070.1000RS
	75.9	Blue	1000 µL regular tips, bevel point, autoclavable		1000 / bag	10 bags	3190.1000B
		Blue	1000 µL regular tips, bevel point, autoclavable		6x 96 / rack	8 packs	3190.1000R
		Blue	1000 µL regular tips, bevel point, autoclavable	✓	6x 96 / rack	8 packs	3190.1000RS
	75.9	Natural	1000 µL regular tips, bevel point, graduated, autoclavable		1000 / bag	10 bags	3090.1000B
		Natural	1000 µL regular tips, bevel point, graduated, autoclavable		6x 96 / rack	8 packs	3090.1000R
	72.0	Blue	1000 µL regular tips, autoclavable		4x 250 / bag	N/A	319.1000B
		Blue	1000 µL regular tips, autoclavable		8x 60 / rack	N/A	319.1000R

1000 µL Filter Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	88.9	Natural	1000 µL universal tips, fine pointed	✓	6x 96 / rack	8 packs	3070.1000FRS

1000 µL Low Protein Binding Tips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	88.9	Natural	1000 µL universal tips, low binding, fine pointed, autoclavable		6x 96 / rack	8 packs	3070.1000PR
		Natural	10 µL ultra-microtips, fine pointed, autoclavable	✓	6x 96 / rack	8 packs	3070.1000PRS

2, 5, 10 mL Macrotips	Size mm	Description			Packaging		Cat. No.
		Colour	Style	Sterile	Pack	Case	
	115.5	Natural	2 mL macrotips, autoclavable		2x 250 / bag	N/A	312.02
	123.2	Natural	5 mL macrotips, graduated, autoclavable		2x 250 / bag	N/A	312.05B
		Natural	5 mL macrotips, graduated, autoclavable		2x 50 / rack	N/A	312.05R
	150.0	Natural	10 mL macrotips, autoclavable		3x 100 / bag	N/A	312.10
		Natural	Empty rack for 24x 10 ml macrotips, autoclavable		4 / pack	N/A	312.10ER

B = bag, E = empty, F = filter, G = gel, L = long, P = low binding, R = rack, S = sterile, T = *tipfill™* refill system

Socorex® compatibility chart

Chart indicates optimal combination of Qualitix® and Socorex® pipettes.

Chart indicates optimal combination of Qualitix® and Socorex® pipettes.			Acura® manual and Acura® electro																						
			815					825 / 826 / 926							835 / 936			855 / 956							
			Up to 10 µL	10 µL / 15 / 20 µL	25 to 100 µL	120 to 200 µL	250 to 1000 µL	0.1 - 2 µL	0.5 - 10 µL	1 - 10 µL (Y)	2 - 20 µL	5 - 50 µL	10 - 100 µL	20 - 200 µL	100 - 1000 µL	0.2 - 2 mL	0.5 - 5 mL	1- 10 mL	0.5 - 10 µL	5 - 50 µL	10 - 100 µL	20 - 200 µL	40 - 350 µL		
10 µL Microtips	    	Tip ref.																							
		3090	●					●	●									●							
		3020	●					●	●																
		302	● ²						● ²																
		3090	●					●	●										●						
3020	●					●	●																		
20 - 100 µL Microtips			Tip ref.																						
 			3080		●					●	●	● ¹													
			3080			●						●	●	● ¹											
200 µL Microtips			Tip ref.																						
     			3070		●	●	●			●	●	●	●	●					●	●	●				
			3290		●	●	●			●	●	●	●	●						●	●	●			
			328																						
			3080				●				●	●	●	●											
			3070			●	●							●								●			
			3080				●								●										
300 - 350 µL Microtips			Tip ref.																						
 			308				●					●	●								●				
			308				●							●								●			
1000 µL Microtips			Tip ref.																						
    			3070				●						●												
			3190				●							●											
			3090				●							●											
			319				●							●											
			3070				●								●										
2, 5, 10 mL Macrotips			Tip ref.																						
  			312											●											
			312														●								
			312																●						

¹ Up to max. tip volume ² For volumes > 2 µL

Ordering information - Instruments

Acura® manual 815, fixed volume micropipettes

Volume	Tip style	Cat. No.	Volume	Tip style	Cat. No.	Volume	Tip style	Cat. No.
1 µL	10 µL	815.0001	50 µL	200 µL	815.0050	250 µL	1000 µL	815.0250
5 µL	10 µL	815.0005	60 µL	200 µL	815.0060	300 µL	1000 µL	815.0300
10 µL	10 µL	815.0010	70 µL	200 µL	815.0070	400 µL	1000 µL	815.0400
10Y µL	200 µL	815.0010Y	75 µL	200 µL	815.0075	500 µL	1000 µL	815.0500
15 µL	200 µL	815.0015	80 µL	200 µL	815.0080	600 µL	1000 µL	815.0600
20 µL	200 µL	815.0020	90 µL	200 µL	815.0090	700 µL	1000 µL	815.0700
25 µL	200 µL	815.0025	100 µL	200 µL	815.0100	750 µL	1000 µL	815.0750
30 µL	200 µL	815.0030	120 µL	200 µL	815.0120	800 µL	1000 µL	815.0800
32 µL	200 µL	815.0032	150 µL	200 µL	815.0150	900 µL	1000 µL	815.0900
40 µL	200 µL	815.0040	200 µL	200 µL	815.0200	1000 µL	1000 µL	815.1000

Acura® manual 825 and 826 XS, adjustable volume micropipettes

Volume	Division	Tip style	Cat. No.	
0.1 – 2 µL	0.002 µL	10 µL	825.0002	826.0002
0.5 – 10 µL	0.01 µL	10 µL	825.0010	826.0010
1 – 10 µL	0.01 µL	200 µL	825.0010Y	826.0010Y
2 – 20 µL	0.02 µL	200 µL	825.0020	826.0020
5 – 50 µL	0.1 µL	200 µL	825.0050	826.0050
10 – 100 µL	0.1 µL	200 µL	825.0100	826.0100
20 – 200 µL	0.2 µL	200 µL	825.0200	826.0200
100 – 1000 µL	1 µL	1000 µL	825.1000	826.1000

Acura® manual 835F, fixed volume micropipettes

Volume	Tip style	Cat. No.
2 mL	2 mL	835.F02
2.5 mL	5 mL	835.F02.5
5 mL	5 mL	835.F05
10 mL	10 mL	835.F10

Acura® manual 835, adjustable volume micropipettes

Volume	Division	Tip style	Cat. No.
0.2 – 2 mL	0.002 mL	2 mL	835.02
0.2 – 2 mL	0.002 mL	2 mL/Pasteur	835.02PP
0.5 – 5 mL	0.01 mL	5 mL	835.05
0.5 – 5 mL	0.01 mL	5 mL/Pasteur	835.05PP
1 – 10 mL	0.01 mL	10 mL	835.10

Acura® manual 855, multichannel micropipettes

Volume	Division	Tip style	8-channel/Cat. No.	12-channel/Cat. No.
0.5 – 10 µL	0.01 µL	ultra 10 µL	855.08.010	855.12.010
5 – 50 µL	0.1 µL	200 µL	855.08.050	855.12.050
10 – 100 µL	0.1 µL	200 µL	855.08.100	
20 – 200 µL	0.2 µL	200 µL	855.08.200	855.12.200
40 – 350 µL	0.4 µL	350 µL	855.08.350	855.12.350

Ordering information - Pipette multipacks

Triopack Acura® manual 825 and 835

Pack	Contents (Micro-/Macropipettes)			Cat. No.
C	2 µL (0.1 – 2)	10 µL (0.5 – 10)	50 µL (5 – 50)	825.TRIO.C
J	2 µL (0.1 – 2)	20 µL (2 – 20)	200 µL (20 – 200)	825.TRIO.J
M	10Y µL (1 – 10)	50 µL (5 – 50)	200 µL (20 – 200)	825.TRIO.M
O	10 µL (0.5 – 10)	100 µL (10 – 100)	1000 µL (100 – 1000)	825.TRIO.O
P	10Y µL (1 – 10)	100 µL (10 – 100)	1000 µL (100 – 1000)	825.TRIO.P
S	20 µL (2 – 20)	200 µL (20 – 200)	1000 µL (100 – 1000)	825.TRIO.S
T	50 µL (5 – 50)	200 µL (20 – 200)	1000 µL (100 – 1000)	825.TRIO.T
V	100 µL (10 – 100)	1000 µL (100 – 1000)	5 mL (0.5 – 5)	825.TRIO.V
Y	1000 µL (100 – 1000)	5 mL (0.5 – 5)	10 mL (1 – 10)	835.TRIO.Y

TwixS Pack Acura® manual 826 XS

Pack	Contents (Micropipettes)		Cat. No.
E	2 µL (0.1 – 2)	20 µL (2 – 20)	826.TWX.E
G	10 µL (0.5 – 10)	100 µL (10 – 100)	826.TWX.G
H	20 µL (2 – 20)	200 µL (20 – 200)	826.TWX.H
K	50 µL (5 – 50)	200 µL (20 – 200)	826.TWX.K
N	100 µL (10 – 100)	1000 µL (100 – 1000)	826.TWX.N
Q	200 µL (20 – 200)	1000 µL (100 – 1000)	826.TWX.Q



► Instant calibration

All Acura® manual models include the *swift-set* user calibration system - with integrated key and locking mechanism - making tedious procedures a thing of the past. Easy and precise, it is in full compliance with international standards and recommendations.




► Safety seal sticker

Though easily removed when needed, the switch protection seal is autoclaving resistant. Replacement service stickers (5 ea.) included with each pipette.

Calibration seal stickers

Replacement seals	25 / pk	1.825.025
Replacement seals	100 / pk	1.825.100

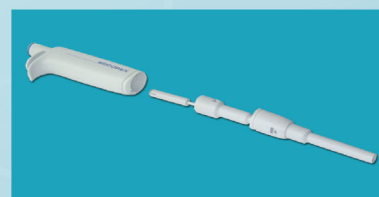

► Minimum maintenance / sterilisation

Easy disassembling for thorough cleaning.

Sterilisation achieved by autoclaving the instrument fully assembled (121°C / 250°F, 20 minutes).



digital reading micropipettes

Ordering information - Accessories

Shelf pipette holder 332

Colour	Packaging	Cat. No.
Transparent blue	1 / pk	320.332
Transparent blue	4 / pk	320.332.4

Twister™ universal 336

Colour	Packaging	Cat. No.
Diamond white	1 / pk	320.336W
Quartz grey	1 / pk	320.336Q
Citrine yellow	1 / pk	320.336Y
Topaz orange	1 / pk	320.336O
Rubin red	1 / pk	320.336R
Emerald green	1 / pk	320.336G
Sapphire blue	1 / pk	320.336B

In-line™ universal 337

Colour	Packaging	Cat. No.
Light grey	1 / pk	320.337G
Pastel rose	1 / pk	320.337R
Ice blue	1 / pk	320.337B
Mint green	1 / pk	320.337M
Vanilla yellow	1 / pk	320.337Y

Static work station 340

Colour	Packaging	Cat. No.
Light grey	1 / pk	320.340

Smartie plunger caps

Colour	Cat. No.	Colour	Cat. No.
Alternative colour selection (6 / pk)			
White	825.710	Rose	825.719
Grey	825.715	Pink	825.720
Vanilla	825.716	Purple	825.721
Lemon	825.712	Ice blue	825.722
Yellow	825.717	Blue	825.713
Orange	825.718	Mint	825.723
Red	825.714	Green	825.711
Smartie mix, assorted colours (14 / pk)			1.825.700

Nozzle protection filter

Description	Packaging	Cat. No.
For 2 and 5 mL models	250 / pk	322.05
For 10 mL model	100 / pk	322.10

Pasteur pipette adapters

Description	Packaging	Cat. No.
For 2 mL model	1 / pk	1.835.631
For 10 mL model	1 / pk	1.835.633

All-purpose reagent reservoirs

Material	Volume	Style	Lid	Sterilized	Packaging	Cat. No.
PP	70 mL	V-shaped	No	Autoclavable	20 / pack	330.01
PP	125 mL	V-shaped	Hinged lid	Autoclavable	20 / pack	330.01L
PS	48 mL	8-channel	No	No	30 / pack	330.08
PS	48 mL	8-channel	No	No	100 / pack	330.08XL
PS	48 mL	8-channel	No	Sterilized	10x 1 / pack	330.08.9
PS	48 mL	12-channel	No	No	30 / pack	330.12
PS	48 mL	12-channel	No	No	100 / pack	330.12XL
PS	48 mL	12-channel	No	Sterilized	10x 1 / pack	330.12.9

Your local distributor

Standards and compliance

All Acura[®] manual models are designed to perform in accordance with national and international standards such as ISO 8655, GMP, GLP, NCCLS.

QC and warranty

Performance values are obtained with bi-dist. water at constant temperature ($\pm 0.5^{\circ}\text{C}$) between 20 and 25°C, in accordance with ISO 8655. Each micropipette bears its own serial number and passes strict performance control attested by an individual QC certificate.

Three-year warranty.



Socorex Isba SA
 Chemin de Champ-Colomb 7a
 1024 Ecublens/Lausanne
 Switzerland
 socorex@socorex.com
 www.socorex.com
 Tel. +41 (0)21 651 6000

B.ACUX.E – C1218